

Test tubes containing phosphate-buffered saline solution



The test tubes containing phosphate-buffered saline solution (PBS) are used as transport containers for swab samples in the veterinary field. Samples can be collected using a swab and absorbed into the solution in the tube.

PBS solution is an isotonic, pH-buffered saline solution designed to maintain the physiological integrity of the sample during transport. PBS solution is not a culture medium.


Functions of the PBS solution:

- ✓ Prevents the smear from drying out,
- ✓ Keeps cells and microorganisms osmotically stable and
- ✓ Has as little chemical effect on the sample as possible

On request, the tubes are available in various sizes and designs. After manufacturing, the products undergo a radiation treatment process. On request, the caps are also available in other colors and, if desired, with individual (barcode) labels.

Item no.	Description	Packing unit
0961 327X	8 ml and 2 ml test tubes containing PBS solution, with individual labels or barcode labels (with tear-off)	On request

Upon request, we can also supply you with compatible swabs, as well as tubes with screw caps designed to hold the swab - which has been broken beforehand and inserted into the tube. During automated processing, the swab is removed along with the tube as soon as it is opened, making the sample material easily accessible for processing. For sample storage, we offer suitable screw caps for resealing.

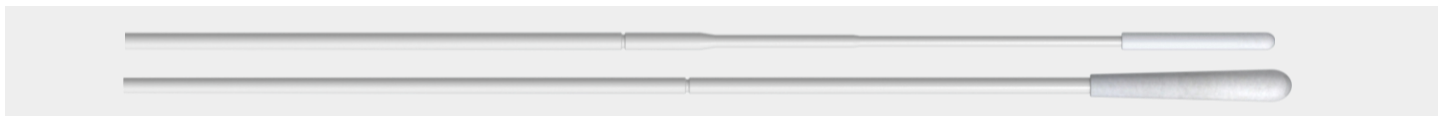

Item no.	Description	Packing unit
056504	Nylon flock fiber swab, nose breakpoint 85 mm, tip-Ø 2.5 mm	On request
056509	Nylon flock fiber swab, throat breakpoint 80 mm, tip-Ø 5 mm	On request